



For an **immediate**
and **natural result**

CYTOCARE
&
REVITACARE Bio-Revitalisation

Reveal the glow and the beauty of your skin

LABORATOIRE
REVITACARE
FRANCE

MADE IN FRANCE



Beauty starts with the quality of the skin

Years go by, stress, hormonal imbalances and external aggressions such as sun, cold and pollution play a negative role on the quality of the skin.

Fine lines and wrinkles appears, skin becomes thinner, duller and loses its elasticity due to the loss of hyaluronic acid naturally present in the skin.

REVITACARE developed a complete range of **SKIN REJUVENATING BOOSTER** based on non cross-linked Hyaluronic Acid to preserve and restore the physiological quality of the skin by promoting the conservation of its elasticity and protecting against external aggressions.





THE SKIN REJUVENATING BOOSTER

by REVITACARE





The complete anti-aging range for an immediate and extremely natural result

The **SKIN REJUVENATING BOOSTER** by REVITACARE offers the most effective solution for all ages to improve in depth the quality of the skin and in surface to give radiance to the skin. This synergistic action promotes the rejuvenation and radiance of the skin for an immediate and natural result.





AN EXCLUSIVE FORMULA BASED ON HYALURONIC ACID AND THE REJUVENATING COMPLEX CT50

-  Fills fine lines,
-  Smooths epidermis,
-  Boosts radiance of skin,
-  Protects from oxidizing agents.

HYALURONIC ACID (2mg - 16mg - 32mg)
REJUVENATING COMPLEX CT50*

In depth action
Preservation of skin
elasticity

FIBROBLASTS

+

On-surface action
Protection against attacks

KERATINOCYTES

=

**DENSIFICATION
OF THE SKIN**

3

different concentrations
for an **optimal efficiency**
at every age and
for every skin type

PREVENTION



CYTOCARE **502**

- ✿ Preserves the youthfulness of the skin
- ✿ Delays the appearance of fine lines
- ✿ Maintains hydration
- ✿ Enhances radiance

CORRECTION



CYTOCARE **516**

- ✿ Smooths early fine lines
- ✿ Increases tonicity and elasticity
- ✿ Restores hydration
- ✿ Improves radiance

REDUCTION




CYTOCARE **532**

- ✿ Corrects skin aging process
 - ✿ Reduces fine lines and wrinkles
 - ✿ Restores hydration, tonicity and elasticity of the skin
 - ✿ Improves radiance
-

* Ammonium Molybdate, Ammonium Metavanadate, Calcium Chloride, Iron Sulfate, Potassium Chloride, Copper Sulfate, Magnesium Chloride, Manganese Sulfate, Sodium Acetate, Sodium Hydrogen Carbonate, Sodium Chloride, Sodium Hydrogen Phosphate, Sodium Metasilicate, Sodium Selenite, Tin Chloride, Zinc Sulfate, Alanine, Arginine, Asparagine, Aspartic Acid, Cysteine, Glutamine, Glutamic Acid, Glycine, Histidine, Isoleucine, Leucine, Lysine, Methionine, Phenylalanine, Proline, Serine, Threonine, Tryptophan, Tyrosine, Valine, Adenine, Biotin, Calcium Pantothenate, Coline Chloride, Folic Acid, Inositol, Nicotinamide, Pyridoxine, Riboflavin, Thiamine, Cyanocobalamin, Deoxythymidine, Glucose, Putrescine, Sodium Pyruvate, Lipoic Acid.

OPTIMIZATION OF ANTIOXIDANT- PROTECTIVE POWER OF HYALURONIC ACID BY CT50

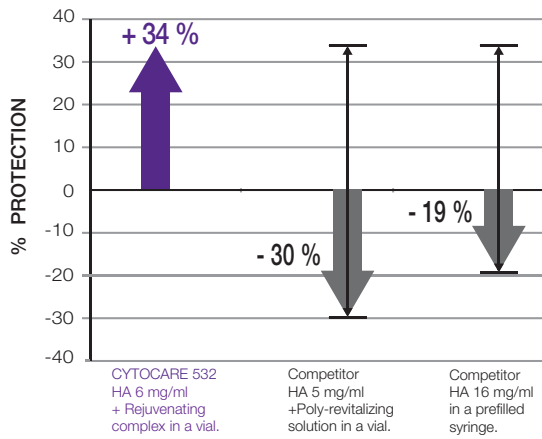


➡ Best antioxidant-protective power *in vivo*.

TABLE 1

Under the experimental conditions of this study, only CYTOCARE has efficiently protected keratinocytes against irradiation, and then is the only one to have antioxidant-protective action.

OPTIMIZATION OF HYDRATING POWER OF HYALURONIC ACID BY CT50

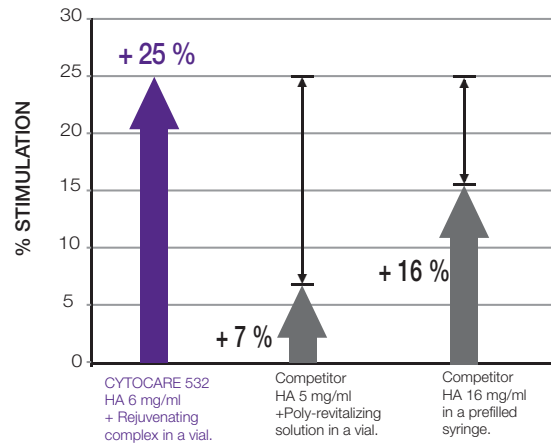

➡ Best hydrating capacity *in vivo*.

TABLE 2

Under the experimental conditions of this study, only CYTOCARE has significantly stimulated sulfated GAGs neosynthesis, which means greater hydrating action.

* *in vitro* studies performed by a French independent laboratory

AN EXCLUSIVE FORMULA BASED ON HYALURONIC ACID AND THE REJUVENATING COMPLEX CT50

VISIBLE RESULTS





CHARACTERISTICS

- ✿ **Box:** 5 vials x 5ml - 10 vials x 5ml,
 - ✿ **Sterilization:** aseptic filtration
 - ✿ **Pyrogen free**
 - ✿ **Ready to use**
 - ✿ **Conservation condition:** 2°C to 30°C
 - ✿ **Close to neutral pH**
 - ✿ **Iso-osmolar solution**
 - ✿ **Injection methods:** manual or gun
 - ✿ **Recommended Protocol:** 3 sessions with 3 weeks intervals, then maintenance of the results by 1 session every 3 months. Protocol and concentration can be adapted depending on the maturity of the skin.
 - ✿ **Injection area:** middle dermis or superficial dermis of the face or the neck.
-

CE 0499



ANTI-AGING TREATMENT FOR MATURE SKIN WITH A NATURAL AND IMMEDIATE RESULT

- 🌸 Fills cutaneous depressions,
- 🌸 Hydrates and Smoothes epidermis,
- 🌸 Lights up the skin.

REVITACARE Bio-Revitalisation, an exclusive combination of high viscosity non cross-linked Hyaluronic Acid and Multivitamins, which provides a spectacular brightness to the skin.

REVITACARE Bio-Revitalisation is particularly recommended to fill crow's feet and perioral folds which are difficult to be treated with a cross-linked Hyaluronic Acid.

COMBINED ACTION OF

Hyaluronic Acid
40 mg

+

Multivitamins*

=

GLOWING SKIN TONE
SMOOTHING OF EPIDERMIS
REDENSIFICATION OF DERMIS

* Ascorbic Acid, Calcium Pantothenate, Ergocalciferol, Nicotinamide, Pyridoxine, Retinol, Riboflavin, Thiamine, Tocopherol, Solutol HS 15.



GLOBAL SATISFACTION*


- **Satisfaction level of patients:**

100% are satisfied
with 65.5% highly satisfied

- **Satisfaction level of practitioners:**

99.5% are satisfied
with 49.5% highly satisfied


- **EVALUATION CRITERIA:**
 - Filling fine lines and wrinkles
 - Hydration
 - Radiance

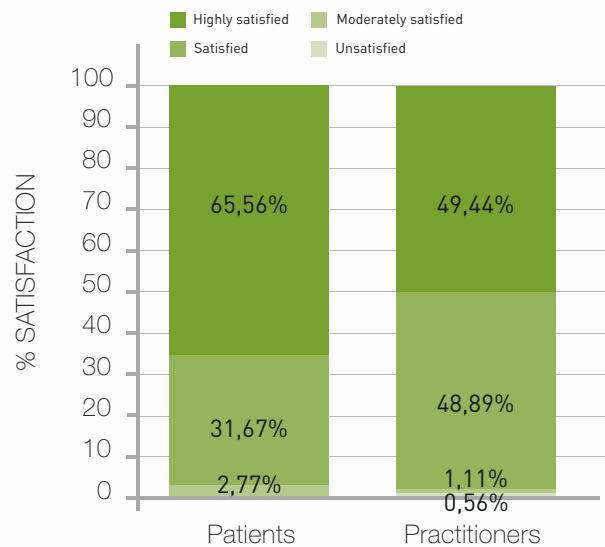

* Global satisfaction level of patients and practitioners at the end of the protocol.
Average duration of follow-up: 13 weeks.
(Retrospective study of 60 patients / 2 centers PROT/CLIN/RV002V01)

CHARACTERISTICS

- ✿ **Box of 2 vials:** Hyaluronic Acid 4ml + Multivitamins 10ml
- ✿ **Sterilization:**
 - Hyaluronic Acid: moist heat
 - Multivitamins: aseptic filtration
- ✿ **Pyrogen free**
- ✿ **Extemporaneously mix inside the syringe** the complete contents of the two vials
- ✿ **Conservation condition:** 2°C to 30°C
- ✿ **Iso-osmolar solution**
- ✿ **Injection methods:** manual or gun
- ✿ **Recommended Protocol:** 3 sessions with 3 weeks intervals, then maintenance of the results by 1 session every 3 months. Protocol and concentration can be adapted depending on the maturity of the skin.
- ✿ **Injection area:** middle dermis or superficial dermis of the face, the neck and the décolleté.

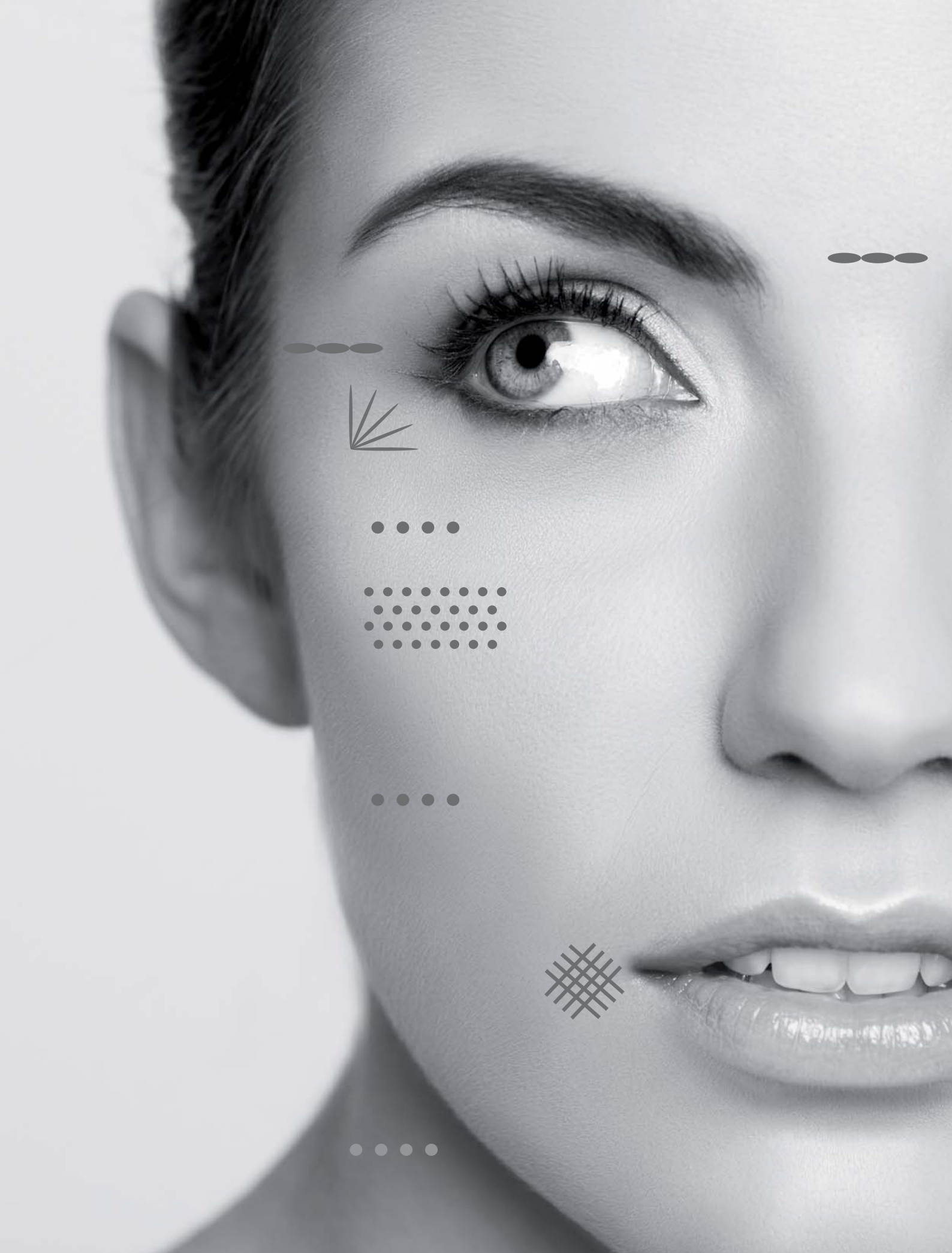

INJECTION TECHNIQUES

CYTOCARE
502

Nappage - point by point - transdermal systems

CYTOCARE
516

Point by point - retracing - fan injection technique - cross-hatching technique - transdermal systems

CYTOCARE
532

Point by point - retracing - fan injection technique - cross-hatching technique - transdermal systems

REVITACARE
Bio-Revitalisation

Nappage - point by point - retracing - fan injection technique - cross-hatching technique

READ CAREFULLY THE LEAFLET BEFORE USING THE PRODUCTS

Needle 13mm/30 G



Fan injection technique



Cross-hatching technique



Retracing

Needle 4mm/30 G



Nappage



Point-by-point





DISTRIBUTOR

REVITACARE is a French laboratory, founded in 2003, specialized in research design and manufacture of implants and products dedicated to cosmetology and aesthetic medicine used in the prevention and correction of the skin aging effects, and the treatment of hair problems.

REVITACARE's products comply with the essential requirements of the Council Directive 93/42/EEC concerning medical devices and its amendments.

QUALITY
SECURITY
EFFICIENCY

- REVITACARE Laboratory: EN ISO 13485 certified.
- Class III Medical Devices.
- HA Quality: biotechnological origin, non cross-linked.
- Microbiological safety.
- Biocompatible in compliance with ISO 10993 international standards: non cytotoxic - non irritant - non sensitizing - non mutagenic.
- Proved efficiency thanks to *in vitro* studies confirmed by clinical data.
- Manufacturing process according to GMP.

LABORATOIRE
REVITACARE
FRANCE

Heliparc - 6 Rue Saint Simon - Parc d'Activités du Vert Galant - 95310 Saint Ouen l'Aumône - FRANCE

Tel: +33 (0) 1 30 37 31 45 - Fax: +33 (0) 1 30 37 88 93 - Email: contact@revitacare.fr

www.revitacare.net